Smart Subscription for CT Scanner Installation Manual

OPERATING DOCUMENTATION

5826427-1EN
Revision 1
© 2019 General Electric Company
All rights reserved.

LEGAL NOTES

TRADEMARKS

Adobe, the Adobe logo, Acrobat, the Acrobat logo, Exchange, and PostScript are trademarks of Adobe Systems Incorporated or its subsidiaries and may be registered in certain jurisdictions.

Microsoft is a registered trademark and Windows is a trademark of Microsoft Corporation.

All other products and their name brands are trademarks of their respective holders.

COPYRIGHTS

All Material Copyrighted (c) 2019 by the GE Company, All rights reserved.

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Damage in Transportation

All packages should be closely examined at time of delivery. If damage is apparent, have notation "damage in shipment" written on all copies of the freight or express bill before delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage MUST be reported to the carrier immediately upon discovery, or in any event, within 14 days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this 14 day period.

To file a report:

- Call 1-800-548-3366 and use option 6.
- Fill out the GIQ workflow for any items missing, damaged, OBF/FOI for in process installs: http://supportcentral.ge.com/ProcessMaps/form_new_request.asp?prod_id=268679&form_id=573167&node_id=1916016&map_id=&reference_id=&reference_type
- Contact your local service coordinator for more information on this process.

Rev. Nov. 10, 2017

Certified Electrical Contractor Statement

All electrical Installations that are preliminary to positioning of the equipment at the site prepared for the equipment shall be performed by licensed electrical contractors. In addition, electrical feeds into the Power Distribution Unit shall be performed by licensed electrical contractors. Other connections between pieces of electrical equipment, calibrations and testing shall be performed by qualified GE Medical personnel. The products involved (and the accompanying electrical installations) are highly sophisticated, and special engineering competence is required. In performing all electrical work on these products, GE will use its own specially trained field engineers. All of GE's electrical work on these products will comply with the requirements of the applicable electrical codes.

The purchaser of GE equipment shall only utilize qualified personnel (i.e., GE's field engineers, personnel of third-party service companies with equivalent training, or licensed electricians) to perform electrical servicing on the equipment.

Lithium Battery Cautionary Statements

RISK OF EXPLOSION.

Danger of explosion if battery is incorrectly replaced.

Replace only with the same or equivalent type recommended by the manufacturer. Discard used batteries according to the manufacturer's instructions.

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Omissions & Errors

Customers, please contact your GE Sales or Service representatives.

GE personnel, please use the Healthcare PQR Process to report all omissions, errors, and defects in this publication.

Contents

1 Readme	1
1.1 Revision History	1
2 Prevace Publication Conventions	3
2.1 Preface	3
2.2 Safety & Hazard Information	3
2.2.1 Text and Character Representation	3
2.2.2 Graphical Representation	4
2.3 Publication Conventions	5
2.3.1 General Paragraph and Character Styles	5
2.3.2 Page Layout	6
2.3.3 Computer Screen Output/Input Text Character Styles	6
2.3.4 Buttons, Switches and Keyboard Inputs (Hard & Soft Keys)	7
3 Procedure	9
3.1 Smart Subscription Installation Flow Chart	9
3.2 Prerequisite	10
3.3 Reference Document	10
3.4 Preparation	10
3.5 Unpacking the ML350 (EES)	10
3.6 Cabling Between ML350 and Switch (EES)	11
3.7 Connect FE Laptop to EES	11
3.8 ML350 Hardware Status Check (EES)	11
3.9 Site Specific Configuration (EES)	11
3.10 Insite (RSvP) Connection (EES)	11
3.11 Check Private IP address of EES (EES)	12
3.12 CT Option Key Install [SmartSubscription -connection] (CT)	12
3.13 Remote Applications Installation (EES)	12
3.14 ML350 Platform Software Status Check (EES)	12
3.15 Cleanup (EES)	13

Page intentionally left blank

1 Readme

1.1 Revision History

Revision	Date	Reason for change
1	July 08, 2019	Initial Release

Page intentionally left blank

2 Prevace Publication Conventions

2.1 Preface

Purpose: Standardize conventions for representing information is a uniform way of communicating information to a reader in a consistent manner. Conventions are used so that the reader can easily recognize the actions or decisions that must be made. There are a number of character and paragraph styles used in this publication to accomplish this task. Please become familiar with them before proceeding forward.

2.2 Safety & Hazard Information

2.2.1 Text and Character Representation

Within this publication, different paragraph and character styles have been used to indicate potential hazards. Paragraph prefixes, such as hazard, caution, danger and warning, are used to identify important safety information. Text (Hazard) styles are applied to the paragraph contents that is applicable to each specific safety statement. Words describe the type of potential hazard that may be encountered and are placed immediately before the paragraph it modifies. Safety information will normally include:

- Type of potential Hazard
- Nature of potential injury
- Causative condition
- How to avoid or correct the causative condition

EXAMPLES OF HAZARD STATEMENTS USED

It's important that you read and understand hazard statements, and not just ignore them.

SEVERE PERSONAL INJURY

DANGER IS USED WHEN A HAZARD EXISTS WHICH WILL CAUSE SEVERE PERSONAL INJURY OR DEATH IF INSTRUCTIONS ARE IGNORED. THEY CAN INCLUDE:

- ELECTROCUTION
- CRUSHING
- RADIATION

SEVERE PERSONAL INJURY

WARNING IS USED WHEN A HAZARD EXISTS WHICH COULD OR CAN CAUSE SERIOUS PERSONAL INJURY OR DEATH IF INSTRUCTIONS ARE IGNORED. THEY CAN INCLUDE:

- POTENTIAL FOR SHOCK
- EXPOSED WIRES
- FAILURE TO TAG AND LOCKOUT SYSTEM POWER COULD ALLOW FOR UN-COMMANDED MOTION.

MINOR PERSONAL INJURY

Caution is used when a hazard exists which can or could cause minor injury to self or others if instructions are ignored. They include for example:

- Loss of critical patient data
- Crush or pinch points
- Sharp objects

PROPERTY DAMAGE ONLY

Notice is used when a hazard is present that can cause property damage but has absolutely no personal injury risk. They can include:

- Disk drive will be destroyed
- Internal mechanical damage, such as to the x-ray tube
- Coasting the rotor through resonance.

2.2.2 Graphical Representation

Important information will always be preceded by the exclamation point

contained within a

triangle, as seen throughout the Preface chapter. In addition to text, several different graphical icons (symbols) may be used to make you aware of specific types of hazards that could possibly cause harm.

2.3 Publication Conventions

2.3.1 *General Paragraph and Character Styles*

Prefixes are used to highlight important non-safety related information. Paragraph prefixes (such as Purpose, Example, Comment or Note) are used to identify important but non-safety related information. Text styles are also applied to text within each paragraph modified by the specific prefix.

EXAMPLES OF PREFIXES USED FOR GENERAL INFORMATION:

Purpose:

Introduces and provides meaning as to the information contained within the chapter, section or subsection (such as used at the beginning this chapter, for example).

NOTE

Conveys information that should be considered important to the reader.

Example:

Used to make the reader aware that the paragraph(s) that follow are examples of information possibly stated previously.


Comment:

Represents “additional” information that may or may not be relevant to your situation.

2.3.2 Page Layout

Publication Part Number & Revision Number

Publication Title

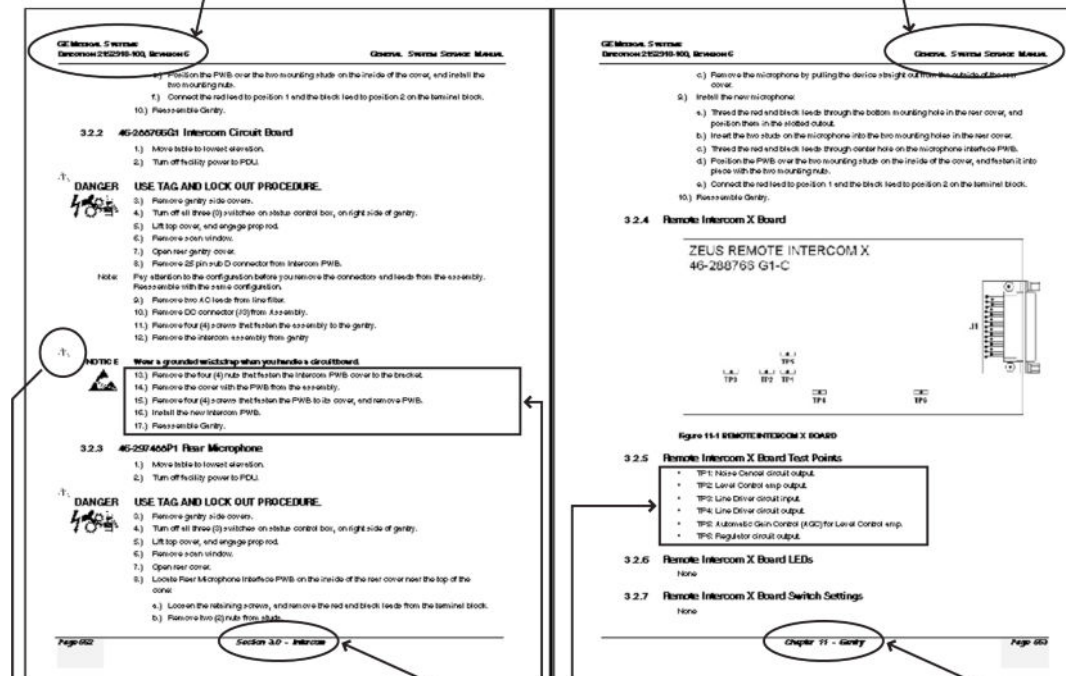The current section and its title are always shown in the footer of the left (even) page.

An exclamation point in a triangle is used to indicate important information to the user.

Paragraphs preceded by **Alphanumeric** characters (e.g. numbers) contain information that must be followed in a **specific order**.

The current chapter and its title are always shown in the footer of the right (odd) page.

Paragraphs preceded by a **symbol** (e.g. bullets) contain information that has **no specific order**.

Headers and footers in this publication are designed to allow you to quickly identify your location. The document part number and revision number appear in every header on every page. Odd numbered page footers indicate the current chapter, its title and current page number. Even numbered page footers show the current section and its title, as well as the current page number.

2.3.3 Computer Screen Output/Input Text Character Styles

Within this publication, mono-spaced character styles (fonts) are used to indicate computer text that is either screen input or output. Mono-spaced fonts, such as `courier`, are used to indicate text direction. When you type at your keyboard, you are generating computer input. Occasionally you will see the math operator “greater-than” and “less-than” symbols used to indicate the start and finish of variable output. When reading text generated by the computer, you are reading it as computer generated output. In addition to direction, characters are italicized (e.g. *italics*) to indicate information specific to your system or site.

Example: Fixed Output

This paragraph's font represents computer generated screen "fixed" output. Its output is fixed from the sense that it does not vary from application to application. It is the most commonly used style used to indicate filenames, paths and text that do not change from system to system. The character style used is a fixed width such as courier.

Example: Variable Output

This paragraph's font represents computer screen output that is "variable". It is used to represent output that varies from application to application or system to system. Variable output is sometimes found placed between greater-than and less-than operators for clarification. For example: <variable_ouput> or <3.45.120.3>. In both cases, the < and > operators are not part of the actual input.

Example: Fixed Input

This paragraph's font represents fixed input. It is computer input that is typed-in via the keyboard. Typed input that does not vary from application to application or system to system. Fixed text the user is required to supply as input. For example: cd /usr/3p

Example: Variable Input

This paragraph's font represents computer input that can vary from application to application or system to system. With variable text, the user is required to supply system dependent input or information. Variable input sometimes is placed between greater-than and less-than operators. For example: <variable_input>. In these cases, the (<>) operators would be dropped prior to input. For example: ypcat hosts | grep <3.45.120.3> would be typed into the computer as

ypcat hosts | grep 3.45.120.3

without the greater-than and less-than operators.

2.3.4 Buttons, Switches and Keyboard Inputs (Hard & Soft Keys)

Different character styles are used to indicate actions requiring the reader to press either a hard or soft button, switch or key. Physical hardware, such as buttons and switches, are called hard keys because they are hard wired or mechanical in nature. A keyboard or on/off switch would be a hard key. Software or computer generated buttons are called soft keys because they are software generated. Software driven menu buttons are an example of such keys. Soft and hard keys are represented differently in this publication.

Example: Hard Keys

A power switch **ON/OFF** or a keyboard key like **ENTER** is indicated by applying a character style that uses both over and under-lined bold text. This is a hard key.

Example: Soft Keys


Whereas the computer **MENU** button that you would click with your mouse or touch with your hand uses over and under-lined regular text. This is a soft key.

Page intentionally left blank

3 Procedure

3.1 Smart Subscription Installation Flow Chart

Figure 1 Smart Subscription installation Flow Chart

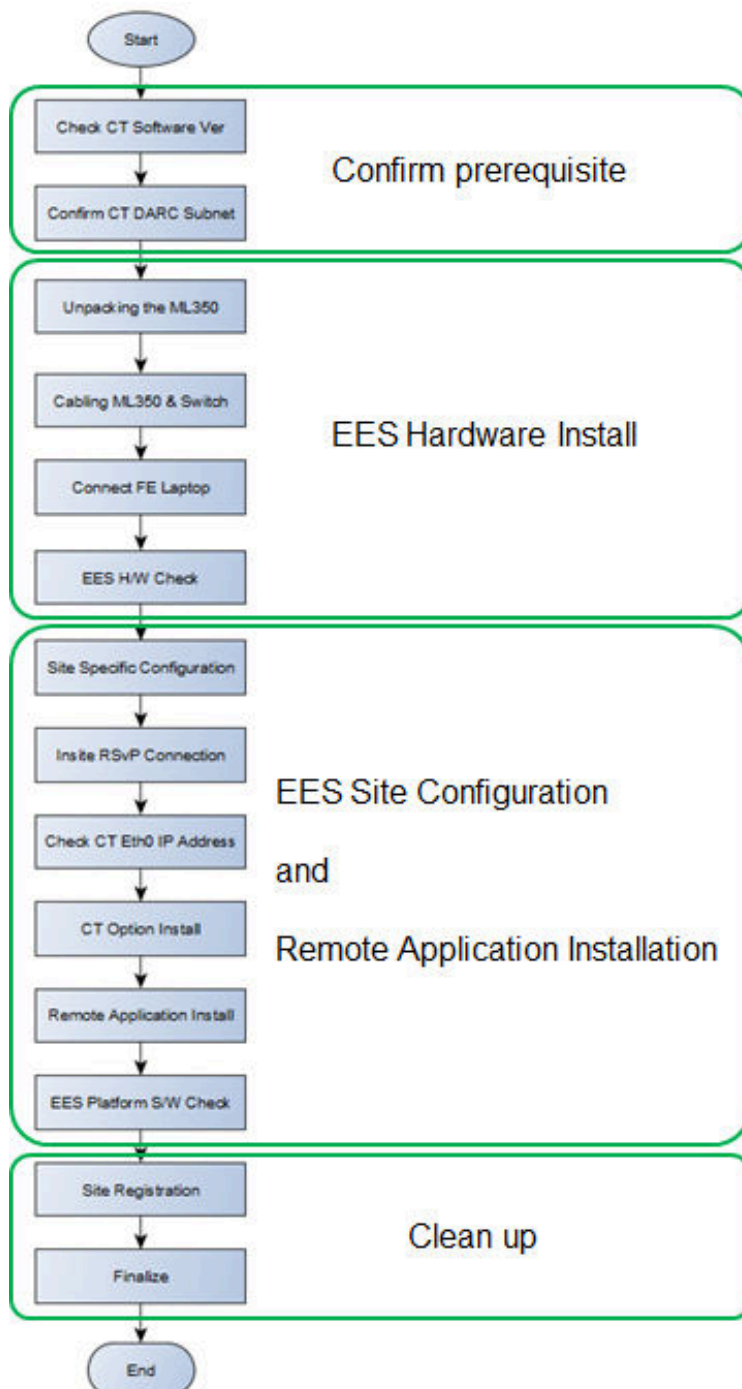3.2 Prerequisite

1. EES Platform Software Installation is completed.
2. 4 consecutive static IP Address are already provided by customer.
3. System ID for EES is required for **RSvP** and **Site** Configuration.
4. One Ethernet cable is required to connect the EES to the FE laptop.

Refer to **[Required Conditions]** of Common Software Procedure in **[Edge Platform Service Guide]**

3.3 Reference Document

- Edge Platform Pre-Installation Manual
- Edge Platform Installation Manual
- Edge Platform Services Guide
- CT Scanner Service Method (Installation → Smart Subscription)

3.4 Preparation

- **Save User Protocols on CT Scanner (CT)**

This step performs User Protocols Backup.

- **Check CT Scanner Software Version (CT)**

This step confirms that the CT scanner software version is compatible with Smart Subscription.

NOTE

If no match, it is necessary to perform the LFC or Install Service Pack before proceeding next step.

- **Confirm Change DARC Subnet is not activate (CT)**

This step is confirmation of the CT Scanner DARC Subnet is not configure **User Defined IP**.

NOTE

If DARC Subnet has been changed already in the Site, it is necessary to EES Platform Software Installation.

Refer to **[Smart Subscription Feature Connection]** in **[Revolution EVO Service Methods]**.

3.5 Unpacking the ML350 (EES)

Two persons are necessary for handling, unpacking and moving the HP ML350G10 server into place.

Refer to **[Unpacking the Hardware Units]** in **[Edge Platform Installation Manual]**.

3.6 Cabling Between ML350 and Switch (EES)

Refer to **[Hardware Installation]** in **[Edge Platform Installation Manual]** for ML350 and Switch.

Refer to **[Smart Subscription Feature Connection]** in **[Revolution EVO Service Methods]** for connect CT Scanner and EES.

3.7 Connect FE Laptop to EES

This step is performing configure the FE Laptop for connect to EES.

Refer to **[Connecting Laptop to Server]** in **[Edge Platform Service Guide]**.

3.8 ML350 Hardware Status Check (EES)

This step performs a check of the ML350 hardware status used by iLO [System Information – Health Summary].

Refer to **[launch iLO Web Interface]** in **[Edge Platform Service Guide]** for login iLO.

- Default Local login name & Password : ctadmin

Refer to **[iLO System Information Selection Screen]** in **[Edge Platform Service Guide]**.

Make sure that the M350 hardware status is no error by using iLO System Information Health Summary.

3.9 Site Specific Configuration (EES)

In this procedure, configure EES according to the Customer environment.

- The setting requires the "sequential IP address" provided by the customer.
- This procedure also includes settings for external gateway, network subnet, DNS, and NTP.

Refer to **[Site Configuration]** in **[Edge Platform installation Manual]**.

NOTE

For **NTP** IP address, be sure to enter the IP address of CT eth0.

* Usually the default is 172.16.0.1

3.10 Insite (RSvP) Connection (EES)

This step is performing the [RSvP] Site Configuration.

Refer to **[RSvP Site Configuration]** in **[Edge Platform Service Guide]**.

NOTE

The [CRM Number] is required in this configuration.

3.11 Check Private IP address of EES (EES)

This step is performing check the two IP address on EES.

Refer to **[Get Server IP]** in **[Edge Platform Installation Manual]**.

Record the following two IP Addresses.

- Edison proxy floating IP (external-net0) _____
- Edison private proxy IP (external-net1) _____

Figure 2 Floating IPs

<input type="checkbox"/>	edison	172.16.0.120	-	external-net1
<input type="checkbox"/>	edison	3.246.33.23	-	external-net0

3.12 CT Option Key Install [SmartSubscription - connection] (CT)

This step performs the [SmartSubscription -connection] CT option installation.

After this option installed, CT scanner and EES will communication will establish.

Refer to **[Smart Subscription Feature Connection]** in **[Revolution EVO Service Methods]**.

3.13 Remote Applications Installation (EES)

The EES is providing some application install procedure.

That is necessary to select the install method by each application.


Refer to each application installation manual for more detail.

Refer to **[Application Software Installation Procedure]** in **[Edge Platform installation Manual]**.

3.14 ML350 Platform Software Status Check (EES)

This step is to confirm the EES platform software is running properly after application install.

Refer to SW T/S part of **[Troubleshooting Flowchart]** in **[Edge Platform Service Guide]**.

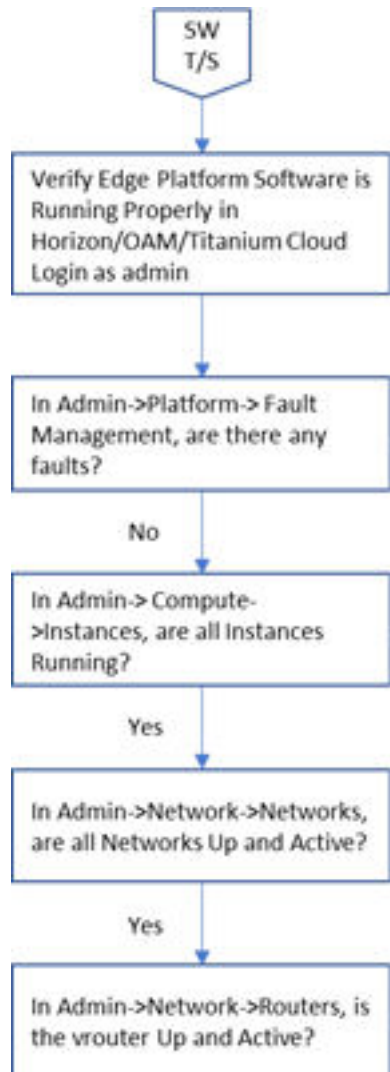Figure 3 Troubleshooting Flowchart

3.15 Cleanup (EES)

Record the site specific parameters using by Site Configuration Sheet.

Refer to [**Appendix X Site Configuration Sheet**] in [**Edge Platform installation Manual**].

Site Registration

Fill out the Site information sheet attached and Send it to GEHCDAAS@he.com.

Figure 4 GEHC Site Installation Information

ITEM	VALUE
Network Details (Site Enterprise IPs)	
iLO	
Site Cidr/subnet	
OAM Public	
Admin Console Public URL	
site gateway	
data-network-range-start IP	
data-network-range-end IP	
Passwords	
wrsroot	
admin	
iLO admin	
Edison-usr	
centos	
BIOS	
Encryption manager	
Applications	
App 1	
App 2	
App 3	
App 4	
App 5	
App 6	
Other Details	
HP Account Username and Password	
Configuration Backup and Restore Password	
Site NAS Access Details	